AL Shqip

Prezentimi

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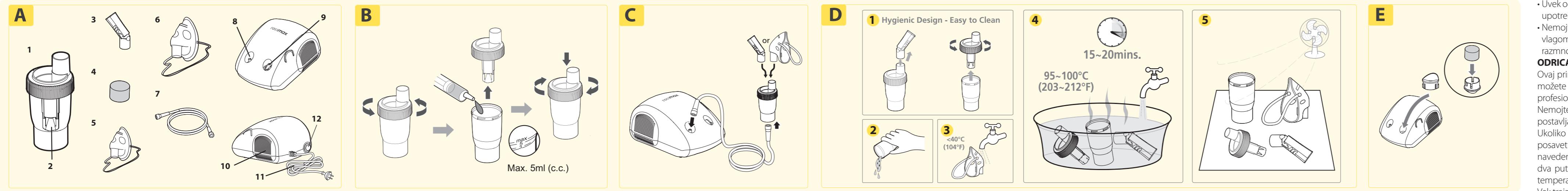
EN English

Introduction

Thank you for purchasing a NE100 Compressor Nebulizer. With proper care and use, your nebulizer will provide you with many years of reliable treatments. This unit operates on standard AC power. Treatments are delivered quickly, safely and conveniently making this unit ideal for all ages in health care structures like hospitals and home care. We encourage you to thoroughly read this guidebook to learn about the features of your nebulizer. Your compressor nebulizer should be used under the supervision of a licensed physician and/or a respiratory therapist. Together with your physician and/or respiratory therapist, you can feel comfortable and confident knowing that you are obtaining the most effective inhalation treatments for your respiratory condition.

NOTE: Your nebulizer is intended for use in treatment of asthma, COPD and other respiratory ailments in which an aerosolized medication is required during therapy. Please consult with your physician and/or pharmacist to determine if your prescription medication is approved for use with this nebulizer. For type, dose and regime of medication follow the instructions of your doctor or respiratory therapist.

This device fulfills the provision of the EC directive 93/42/EEC (Medical Device Directive) and the European Standard EN 13544-1:2007+A1:2009 Respiratory therapy equipment - Part 1: Nebulizing systems and their components.

Please read this manual carefully before use and be sure to keep this manual.

Cautions

Please use general safety precautions when operating your nebulizer. This unit should be used only for its intended purpose as described in this guidebook and with medications only under the supervision and instruction of your physician. Do not use the device in anesthetic or ventilator breathing circuits.

Product cautions

READ THE FOLLOWING BEFORE USING
 • To avoid electrical shock: keep unit away from water.
 • Do not handle the unit of power cord with wet hands.
 • Do not immerse the unit in liquid.
 • Do not use while bathing.
 • Do not reach for a unit that has fallen into water immediately unplug the unit.
 • Do not use the unit if it has any damaged parts (including plug), if it has been submerged in water or dropped. Promptly send the unit for examination and repair.
 • The unit should not be used where flammable gas, oxygen or aerosol spray products are being used.
 • Keep the air vents open. Do not place the unit on a soft surface where the openings can be blocked.
 • If the medication cup is empty, do not attempt to operate the unit.
 • If any abnormality occurs, discontinue to use until the unit has been examined and repaired.
 • The unit should not be left unattended while plugged in.
 • Do not tilt or shake the unit when in operation.
 • Disconnect the unit from the electrical outlet before cleaning, filling and after each use.
 • Do not use attachments unless recommended by the manufacturer.
 • Do not disassemble or attempt to repair the unit.
 • Do not use the device in anesthetic or ventilator breathing circuits.

Operating cautions

• Close adult supervision is highly recommended when the unit is used by children or invalids.
 • Keep your eyes away from the output of medication mist.
 • The maximum capacity of the medication cup is 5 ml and should not be overfilled.
 • Do not use this unit while operating a vehicle.
 • If any discomfort or abnormality occurs, stop using the unit immediately.
 • Do not use the device if the air tube is bent.
 • Pentamidine is not an approved medication for use with this device.
 • Do not use the tubing and/or AC adaptor/power cord for any other purpose than those specified, as they can cause risk of strangulation.
 • Do not service or maintain device while in use.

Storage cautions

• Do not store the unit in direct sunlight, high temperature or humidity.
 • Keep the unit out of reach of small children.
 • Always keep the unit unplugged while not in use.

Cleaning cautions

• Check air filter, nebulizer, mouthpiece and any other optional component before each use. Dirty or worn parts should be replaced.
 • Do not immerse the unit in water. It may damage the unit.
 • Disconnect the unit from the electrical outlet before cleaning.
 • Clean all necessary parts after each use as instructed in this guidebook.
 • Always dispose of any remaining medication in the medication cup after each use. Use fresh medication each time you use the device.
 • Do not store the air tube with moisture or medication remaining in the air tube. This could result in infection as a result of bacteria.

MEDICAL DISCLAIMER:
 This manual and product are not meant to be a substitute for advice provided by your doctor or other medical professionals.

Don't use the information contained herein or this product for diagnosing or treating a health problem or prescribing any medication. If you have or suspect that you have a medical problem, promptly consult your doctor.

Durable periods are as follows, provided the product is used to nebulize 2ml of medication 2 times a day for 8 minutes each time at room temperature (25°C).

Durable period may vary depending on usage environment.

Product specifications

Power	AC 230V/50Hz or AC 220V/60Hz or AC 110V/60Hz
Power Consumption	< 130W
Sound Level	< 60 dB(A (1 meter away from NE100))
Compressor Pressure Range	> 29 psi (200 kPa)
Operating Pressure Range	> 15 psi (103 kPa)
Operating Flow Range	> 3.5 lpm
Operating Temperature Range	10°C to 40°C (50°F to 104°F)
Operating Humidity Range	10 ~ 90% RH
Operating Atmospheric Pressure Range	700-1060 hPa
Storage Humidity Range	10 ~ 90% RH
Storage Atmospheric Pressure Range	700-1060 hPa
Weight	1360g (without accessories)
Particle Size (MMAD)	< 2.6μm
Average Nebulization Rate	> 0.35 ml/min
Standard Accessories	Nebulizer Kit, Air Tube, Mouthpiece, Filters (5pcs), Adult and child masks

* Subject to technical modification without prior notice.

* Performance may vary with drugs such as suspensions or high viscosity. See drug supplier's data sheet for further details.

A Product identification

1. Nebulizer 2. Nozzle 3. Angled Mouthpiece

Kit 4. Air Filter 5. Child Mask 6. Adult Mask

7. Air Tube 8. Air Output 9. Air Filter Slot

10. Air Intake 11. Power Cord 12. Power Switch

B Assembling your nebulizer kit

Follow the cleaning instructions in this guidebook under "Cleaning technique" prior to using your nebulizer for the first time or after it has been stored for an extended period of time.

REMEMBER: Always unplug the compressor and make sure the power-switch is turned to the "OFF" position before cleaning, assembling and before or after each use.

1. Place the compressor on a flat, stable surface within reach.

2. Gently twist and pull straight up on the lid of the nebulizer kit to separate into two parts (medication cup and cover).

3. The stem inside the medication cup inserts into the tube of the nozzle.

4. Add the prescribed amount of medication to the medication cup.

5. Reassemble the nebulizer kit by carefully twisting the medication cup and cover together. Make sure that the two parts fit securely.

WARNING: The symbol on this product means that it's an electronic product and following the European directive 2012/19/EU the electronic products have to be disposed on your local recycling centre for safe treatment.

C Operating your nebulizer

The nebulizer kit is operable at up to a 45° angle. If the angle is greater than 45°, no aerosol will be generated.

1. Attach one end of the tube connector to the air output.

2. Carefully attach the opposite end of the air tubing connector to the stem at the base of the nebulizer kit.

3. Attach the angled mouthpiece or mask to the top of the nebulizer kit.

4. The capacity of the medication cup is 2-5 ml.

NOTE: A 30-minute interval is recommended after each use. The compressor will automatically shut off if it becomes overheated.

If / when this happens, immediately:

1. Press the power-switch to the "OFF" position.

2. Unplug the power cord from the outlet.

3. Allow the motor to cool for 30 minutes.

Before restarting the unit, make sure that the air vents are not obstructed.

D Cleaning procedures

Rinsing technique (performed after each treatment or before first use):

1. Unplug the unit from the power source.

2. Allow the unit to completely cool.

3. Carefully detach the air tubing from the nebulizer and pour out any remaining medication.

4. Follow the cleaning procedures provided in this guidebook.

E Cleaning procedures

Rinsing technique (performed after each treatment or before first use):

1. Make sure that the power-switch has been turned to the "OFF" position and the unit has been disconnected from the power source.

2. Disconnect the air tube from the nebulizer device.

3. Gently twist and pull up the cover of the nebulizer kit to open and separate.

F Cleaning procedures

Rinsing technique (performed after each treatment or before first use):

1. Make sure that the power-switch has been turned to the "OFF" position and the unit has been disconnected from the power source.

2. Disconnect the air tube from the nebulizer device.

3. Gently twist and pull up the cover of the nebulizer kit to open and separate.

G Cleaning procedures

Rinsing technique (performed after each treatment or before first use):

1. Rinse the nebulizer kit and components with hot tap water.

2. Dry with clean towels or completely air dry.

3. Reassemble the nebulizer kit.

4. Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

H Cleaning the air filter

It is important to change the air filter approximately when the air filter turns gray. It is recommended to change air filter every 2 months.

1. Remove the air filter cover by gently pulling forward.

2. Discard the gray filter.

3. Replace with a new, clean air filter.

4. Securely re-attach the air filter cover to the unit.

I Changing the air filter

It is important to change the air filter approximately when the air filter turns gray. It is recommended to change air filter every 2 months.

1. Remove the air filter cover by gently pulling forward.

2. Discard the gray filter.

3. Replace with a new, clean air filter.

4. Securely re-attach the air filter cover to the unit.

J Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

K Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

L Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

M Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

N Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

O Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

P Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

Q Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

R Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

S Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

T Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

U Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

V Cleaning the compressor

Wipe the compressor daily using a soft cloth.

NOTE: Any other form of cleaning or cleaning agents may damage the finish of the unit.

W Cleaning the compressor

Wipe the compressor daily using a soft cloth.